BRAUN

ExactFit® 4

Upper arm blood pressure monitor

BUA6450

Braun ExactFit® 4

BUA6450US

Upper arm blood pressure monitor

Contact Us

Questions or Comments?

Call us toll-free at: 1-800-327-7226

Contact us at: www.BraunHealthcare.com/contact-us

Or visit our website at: www.BraunHealthcare.com

Please be sure to specify the model number.

Guangdong Transtek Medical
Electronics Co., Ltd.

0123

Zone A, No.105, Dongli Road, Torch
Development District,
528437 Zhongshan, Guangdong, China

Certain trademarks used under license from
The Procter & Gamble Company or its affiliates.

ExactFit® is a trademark of Kaz Europe Sàrl.

Distributed by Kaz USA, Inc.

A Helen of Troy Company
400 Donald Lynch Blvd., Suite 300,
Marlborough, MA 01752


Made in China

© 2026 All rights reserved.

MD REF BUA6450

P/N A014502R0

25MAR26

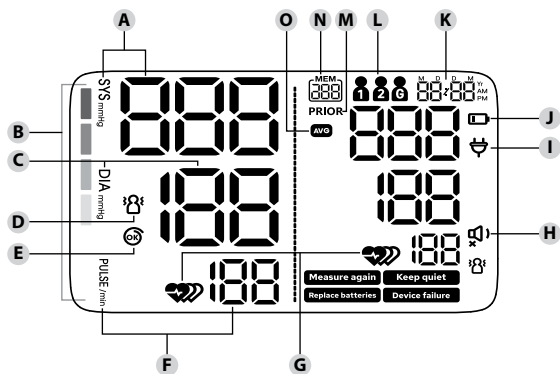Product description

- 1 LCD display
- 2 Colour-coded risk indicator
- 3 Power/Measurement (Start/Stop/Set) button
- 4 Memory/Backward button
- 5 User/Forward button
- 6 Type-C port
- 7 Battery compartment
- 8 Connector
- 9 Arm cuff
- 10 Air hose

Display description

- A Systolic blood pressure reading
- B Colour-coded risk indicator
- C Diastolic blood pressure reading
- D Body motion detection symbol 
Appears when talking, moving, or shaking of the arm with the cuff on is detected during a measurement. NOTE: The measured blood pressure reading may not be accurate when this symbol is displayed with the reading.
- E Cuff fit guide symbol 
Appears when the cuff is wrapped well.
- F Pulse display
Pulse rate appears after the measurement.
- G Pulse rate symbol 
Flashes when detected during a measurement.
- Irregular heartbeat symbol 
Appears when irregular pulse rate is detected during a measurement.
- H Voice symbol 
- I Adapter power mode/Adapter error symbol 
Appears when the power is supplied by AC adapter. Or appears along with "bAt H" to indicate the voltage of the blood pressure monitor is high.
- J Low battery symbol 
Appears when the battery is low along with "bAt Lo".
- K Date/Time display
- L User ID symbol 
The user ID number appears when the monitor is operated by the selected user.
- M Previous measurement symbol **PRIOR**
Indicates the previous measurement.

N Memory symbol

Indicates that the device is in memory mode and displays the number of the stored reading (up to 250).

O Average value symbol

Appears when viewing the average value under the memory mode.

Intended use and contraindications

This Blood Pressure Monitor is intended for use in measuring blood pressure and pulse rate in patients with arm circumferences from 8.6 to 17.7 inch, and is intended for adult population, including those with conditions of diabetes. It is intended for indoor use only.

This blood pressure monitor is not intended for clinical use in a professional environment and is for home usage only. The use of this blood pressure monitor is not intended as a substitute for a consultation with your doctor.

Measurement principle

This product uses the Oscillometric Measuring method to detect blood pressure. Before every measurement, the unit establishes a “zero pressure” equivalent to the atmospheric pressure. Then it starts inflating the arm cuff, meanwhile, the unit detects pressure oscillations generated by beat-to-beat pulsatile, which is used to determine the systolic and diastolic pressure, and also pulse rate.

Receiving and inspecting your monitor

Check that the device packaging has not been tampered with and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

What you should know about blood pressure

Blood pressure constantly changes throughout the day. It rises in the early morning and declines during the late morning. Blood pressure rises again in the afternoon and finally drops to a low level at night. It may vary in a short period of time. Therefore, readings from successive measurements can fluctuate. This device will display your blood pressure reading referring to two values: systolic and diastolic, through the principle of pressure oscillation. The systolic blood pressure (the upper number) indicates how much pressure your blood is exerting against your artery walls when the heart beats. The diastolic blood pressure (the lower number) indicates how much pressure your blood is exerting against your artery walls while the heart is resting between beats. The pulse value is also displayed after each measurement. Blood pressure measured in a doctor's office only provides a momentary value. Repeated measurements at home better reflect one's actual blood pressure values under everyday conditions. Moreover, many people have a different blood pressure when they measure at home, because they tend to be more relaxed than when in the doctor's office. Regular blood pressure measurements taken at home can provide your doctor with valuable information on your normal blood pressure values under «everyday» conditions.

The AHA (American Heart Association) has set up the following standard blood pressure values when measured with a resting pulse at home:

Blood pressure (mmHg)	Normal	Elevated	High Blood Pressure Stage 1	High Blood Pressure Stage 2	Hypertensive Crisis
SYS = systole (upper value)	Less than 120	120-129	130-139	140-180	Higher than 180
DIA = diastole (lower value)	Less than 80	Less than 80	80-89	90-120	Higher than 120

 **CAUTION:** If you are under medical treatment or taking any medication, please consult your doctor first. **DO NOT** adjust medication based on measurement results from this blood pressure monitor. Take medication as prescribed by your doctor. Only a doctor is qualified to diagnose and treat high blood pressure.

Safety information and important precautions

To ensure accurate measuring results, read the complete instructions for use carefully.

WARNING:

If your blood pressure reading is repeatedly in the moderate or severe hypertension range - please consult your doctor or follow your doctor's current recommendation.

DO NOT self-diagnose the measurement results and start treatment by yourself. The measurement results given by this device is not a diagnosis. **ALWAYS** consult your doctor for evaluation of the results and treatment.

DO NOT adjust medication based on readings from this blood pressure monitor. Take medication as prescribed by your physician. **ONLY** a physician is qualified to diagnose and treat high blood pressure.

DO NOT apply the cuff on an arm that has an intravenous drip or a blood transfusion attached.

DO NOT kink, fold, stretch, compress, or otherwise deform the tube during measuring, as the cuff pressure might continuously increase, which could prevent blood flow and result injury.

Measuring blood pressure too frequently can cause injury to the user due to blood flow interference.

DO NOT apply the cuff over a wound as this can cause further injury.

DO NOT apply cuff to areas on patient where skin is delicate or damaged. Check cuff site frequently for irritation.

DO NOT apply the cuff to a limb on which there is intravascular access or on which an arterio-venous (A-V) shunt is present. This can cause temporary interference of blood flow and potential injury to the user.

DO NOT place the cuff on the arm on the same side of a mastectomy (especially when lymph nodes have been removed). it is recommended to take measurements on the unaffected side.

DO NOT use the device while simultaneously using other medical electrical (ME) equipment on the same limb. One or both devices could temporarily stop functioning if you try to use them at the same time.

Please check (for example, by observation of the limb concerned) that the operation of the device does not result in prolonged impairment of patient blood circulation.

On the rare occasion of a fault causing the cuff to remain fully inflated during measurement, loosen and remove the cuff immediately. Prolonged high pressure applied to the arm (cuff pressure >300 mmHg or constant pressure >15 mmHg for more than 3 minutes) might lead to bruising and discolored skin.

DO NOT use this device with high-frequency (HF) surgical equipment at the same time.

This device is not used in conjunction with oxygen rich environments, not intended for use with flammable anaesthetics, not intended for use in conjunction with flammable agents.

Excessive cuff tube lengths could cause strangulation if you don't manage them properly.

DO NOT touch output of the batteries/adaptor and the user simultaneously.

The power cord is considered the disconnect device for isolating this equipment from supply mains.

DO NOT position the equipment so that it is difficult to reach or disconnect.

DO NOT use this device if you are allergic to polyester, nylon, or plastic. Please use the original factory cuff to ensure accuracy.

Using unapproved accessories might cause damage to the unit and injure users.

If you experience discomfort during a measurement, such as pain in the arm or other complaints, press the Power button immediately to release the air from the cuff.

DO NOT use the device while under maintenance, or being serviced.

Sensor degradation or looseness may reduce performance of device or cause other problems.

The air tube poses a risk of strangulation. Furthermore, the small parts of product and batteries present a choking hazard if swallowed. They should therefore always be kept away from infants/children.

CAUTION:

DO NOT attempt to repair the unit yourself if it malfunctions, contact consumer service.

NEVER open the device! This will make the manufacturer's warranty invalid and could adversely affect the device's accuracy!

Store your device, cuff and adaptor in a clean and dry place, protect it against extreme moisture, heat, lint, dust and direct sunlight. Never place any heavy objects on it.

Avoid compression or restriction of the connection tubing (in operation and storing). Continuous cuff pressure due to the connection tubing kinking can affect blood flow and cause possible injury to the user.

PRECAUTION:

This device is portable, but it is not intended for use during patient transport. This device is not suitable for continuous monitoring during medical emergencies or operations.

This device is intended for non-invasive measuring and monitoring of arterial blood pressure. It is not intended for use on extremities other than the arm, or for any purpose other than obtaining a blood pressure measurement.

This device is for adult use, including diabetic patients. **DO NOT** use this device on neonates or infants. **DO NOT** use it on children and adolescents unless otherwise instructed by a medical professional. **DO NOT** use it on patients wearing or implanting electronic devices, such as pacemakers, defibrillators.

Consult with your physician before using this monitor if you suffer from the following conditions: common arrhythmias such as premature ventricular beats or atrial fibrillation; peripheral arterial disease, implantation with cardiac pacemaker; undergoing intravascular therapy; arteriovenous shunt or mastectomy. Please note that any of these conditions may affect measurement readings, in addition to patient motion, trembling or shivering.

If you are taking medication, consult your physician to determine the proper time to measure your blood pressure.

This device may be used only for the intended use described in this manual, the manufacturer shall have no liability for any incidental, consequential, or special damages caused by intentional misuse or abuse.

Please use the device under the environment which is provided in the user manual. Otherwise, the performance and lifetime of the device will be impacted and reduced.

The device may require up to 30 minutes to warm up/cool down from the minimum/maximum storage temperature before it is ready for use.

The blood pressure monitor, its adapter, and the cuff are suitable for use within the patient environment.

DO NOT wash the cuff in a washing machine or dishwasher!

NOTICE:

You can use this device to take your own measurement, no third-party operator is required.

Adapter is specified as a part of ME EQUIPMENT.

When requested, test method will be made available by manufacturer or distributor.


The expected lifetime of the cuff may vary by the frequency of washing, skin condition, and storage state.

Please report to the manufacturer and the competent authority of the Member State in which you are established about any serious incident that has occurred in relation to this device.

Inserting batteries


- Slide off the battery cover.
- Install or replace 4 AA (LR6 1.5V) size batteries according to the polarity indications inside the battery compartment.
- Place back the battery cover.

Any time the battery is low, it will display the icon  and "bAt Lo" and then power off automatically after about 5 seconds.

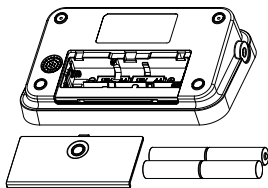Replace the batteries:

- When the low battery symbol appears.
(see the LCD display on the right picture)
- When any button is pressed and nothing is displayed on the screen.

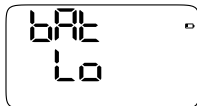⚠ CAUTION:

- New and used batteries, or different types of batteries shall not be used together. Remove batteries if the device is not likely to be used for some time.
- Do not heat or deform the batteries, or dispose of them in fire.

Only discard empty batteries. They should not be disposed of in the household waste, but at appropriate collection sites or at your retailer.


Choice of power supply

- Battery powered mode: 6VDC = 4×AA (LR6 1.5V) batteries OR
- AC adapter powered mode: 5V = 1A

⚠ CAUTION:

In order to maintain performance and protect your device, please use the right batteries and approved AC power adapter which complies with local safety standard.

Please unplug the adapter to disconnect the device from mains power.

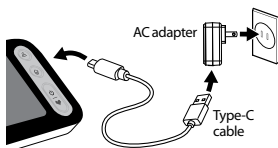Key rules for accurate blood pressure measurement

Always take readings at the same time of day, under the same conditions.

Do not take a measurement within 30 minutes of smoking, drinking coffee or tea or any form of exertion. These factors will influence the measurement.

Ensure at least 5 minutes of quiet rest before measurements; if possible, 10 minutes of rest are recommended. Try to always use the same resting time before taking a measurement.

Always measure on the same arm. The left arm is recommended. Do not move or speak. Sit comfortably, keep your legs uncrossed and feet flat on the floor, back and arm supported.

Place the cuff in direct contact with the skin. Do not wrap the cover over clothing.

If you want to take multiple measurements, wait approximately 5 minutes before repeating.

Applying the arm cuff

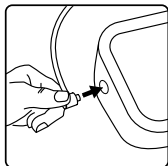
Fig. 1

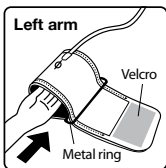
Fig. 2

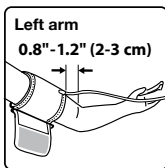
Fig. 3

1. Plug the air hose into the connector (Fig. 1).
2. The cuff is located correctly when the velcro is on the outside of the cuff and the metal ring does not touch the skin (Fig. 2).
3. Put your arm through the cuff loop. The bottom of the cuff should be approximately 0.8-1.2 inches above elbow (Fig. 3). The tube should lie over the brachial artery on the inside of the arm (Fig. 4).

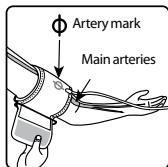
Fig. 4

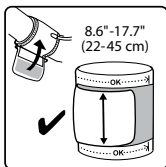
Fig. 5

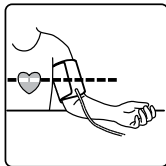
Fig. 6

4. Pull the cuff so that the top and bottom edges are tightened around your arm, make sure not to wrap the cuff very tightly. Allow space for one finger to fit between the cuff and the arm.
5. When the cuff is positioned properly, press the velcro firmly against the cuff.
6. This cuff is ready for use if the arrow heads fall within the range on the other side of the cuff when the cuff is tightened around your arm (Fig. 5).
7. Sit on a chair and place your arm on the table so that the cuff is at the same level as your heart (Fig. 6).

Taking a measurement

1. Sit down with straight back with your feet flat on the floor.
2. Ensure the cuff is wrapped snugly (per the "Applying the arm cuff" section above, Fig. 3-6).
3. Do not move the unit during measurement, or the proper measurement will not be achieved.

4. Press the Power/Measurement button  once to switch on the device. Once powered on, you can navigate the memory and select the desired user. Press the Power/Measurement button  again to start a measurement. While measurement is in progress, do not move or speak.

Note: If necessary to interrupt a measurement, press the Power/Measurement button  at any time. The device immediately cancels the measurement, lowers cuff pressure and automatically turns off.

5. The cuff will first inflate and then deflate. At the end of the measurement the reading is displayed:
Upper (systolic) value
Lower (diastolic) value
Pulse
6. Press the Power/Measurement button  to turn the product off. Otherwise the product will turn off automatically after 30 seconds.
7. Test results will be saved into the memory automatically (see section storing measurement data).

This device is not intended as a substitute for regular check-ups by your doctor, please continue to visit your doctor on a regular basis for a professional reading.

Irregular heartbeat detector

An irregular pulse rate is detected when a pulse rate rhythm varies while the unit is measuring the systolic and diastolic blood pressure. During each measurement, the monitor records all the pulse intervals and calculate the average; if there are two or more pulse intervals, the difference between each interval and the average is more than the average value of $\pm 25\%$, or there are four or more pulse intervals, the difference between each interval and the average is more than the average value of $\pm 15\%$, the irregular heartbeat symbol appears on the display when the measurement results are appeared.

CAUTION:

The appearance of the IHB icon  indicates that a pulse irregularity consistent with an irregular pulse rate was detected during measurement. Usually this is NOT a cause for concern. However, if the symbol appears often, we recommend you seek medical advice. Please note that this feature does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

Color Coded Guidance

The instant read color indicator follows the AHA recommendations for interpreting readings and is described in the table below.

	Green: normal blood pressure		Your blood pressure reading is normal, and values are below 120 mmHg systolic and below 80 mmHg diastolic.
	Yellow: elevated		Your blood pressure reading is in the mild hypertension range with values between 120-129 mmHg systolic and/or less than 80 mmHg diastolic.
	Orange: moderate hypertension		Your blood pressure reading is in the moderate hypertension range with values between 130-139 mmHg systolic and/or 80-89 mmHg diastolic.
	Red: severe hypertension	Warning	Your blood pressure reading is in the severe hypertension range with values of 140 mmHg systolic and above and/or 90 mmHg diastolic and above. If you get readings in this range, consult your doctor straight away.
	Red (Text Flashing): hypertensive crisis	Warning	Your blood pressure reading is in the hypertensive crisis range with values a systolic blood pressure value higher than 180 mmHg and / or a diastolic blood pressure value higher than 120 mmHg. If you find you have one or two readings in this range, consult your doctor right away.

Setting month, date and time

It is important to set the date and time before using your blood pressure monitor for the first time, so that a correct time stamp can be assigned to each record that is stored in the memory.

Note: If using ONLY the power cord, the date and time must be re-set each time the device is unplugged from the power supply. If the power cord remains connected, the date and time will be maintained. If you insert the batteries, the date and time will be stored, you can continue to use the power cord as the main power source.

(The setting range of the year: 2025–2055, Time format: 24H/12H)

- When the device is off, press and hold the Power/Measurement button , and enter the [Year] setting.
Press the User or Memory button to change the year.
- Press the Power/Measurement button  to confirm the [Year], then the date format will flash. Press the User or Memory switch the date format between [month/day] and [day/month].

20 25

20 26


1.

20 25

20 26

2.

- Press the Power/Measurement button  to confirm the date format, then the [Month] will flash. Repeat the same steps to set the [Month] and [Day].
- Repeat the same steps to set the time format between [24H] and [12H].
- Repeat the same steps to set the [Hour] and [Minute].
- After setting the time, the LCD will display all the settings you have done, and then it will display "do nE" and the device will turn off automatically.

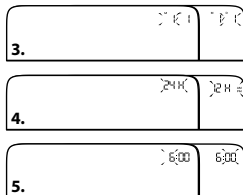
Note: The device can also enter the setup mode automatically when the battery is installed for the first time.

Selecting a user

There are 3 user ID    available. The user  and , each with 250 memory spaces, are designed for 2 different people to save the measured values separately. The user , with no memory space, is reserved for guest.


- Press Power/Measurement button  to turn on device.
- Press the User button to select user 1, user 2, or guest. All readings will be stored for the selected user.
- Press the Power/Measurement button  to start measurement for the selected user. No readings are stored for guest user.

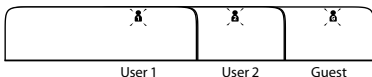Press the Memory button to confirm the selected user ID, it will enter the memory query automatically.

Start a measurement

- Press the Power/Measurement button  to turn on the device. Confirm the correct user is showing on screen.
- Press the Power/Measurement button  again to start the measurement. Remain still and do not talk until the measurement is complete.

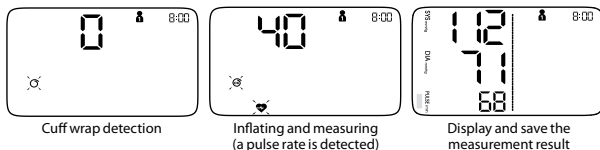Note: Any time, to stop the measurement, press the Power/Measurement button . If the selected user has previous readings stored, it will display its prior measurement after the measurement.

- Press the Power/Measurement button  to turn off the monitor, or it will shut off within 30 seconds.

Note: Both user  and  can store maximum 250 readings. When you pass that limit, every time you take the measurement, the monitor will prompt "FULL" first and the oldest reading will drop from the memory list after the measurement.

- About the irregular pulse rate and excessive body motion during the measurement:

During a measurement, If an irregular pulse rate is detected, the symbol  will display in the measurement result.


During a measurement, when the excessive body motion, especially of the arm the cuff is worn is detected, the symbol  will flash about 5 seconds and check again. If it is no longer detected, the symbol  will disappear; If still detected, the symbol will final display in the measurement result.



Note: The measured blood pressure reading may not be accurate if this symbol is displayed.

Storing measurement data


- Press the Power/Measurement button  to switch the device on, the user ID will start blinking. Press the User button to switch between users. Press the Memory button to select the desired user.
- If the last 3 readings were taken within 30 minutes, the LCD will display the average value of last 3 readings on the right and display the last measurement on the left. If there are fewer than three readings stored in the memory, the latest measurement will be displayed on the right instead.
- Press the User or Memory button to display the next measurement.

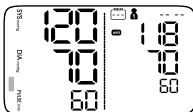User 1 Example

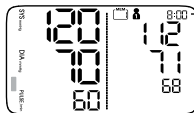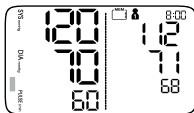The Time and Date will display alternately.

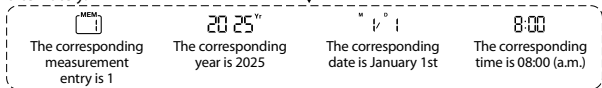Note: The latest measurement (1) is shown first. Each new measurement is assigned to the first (1) measurement. All other readings are pushed back one digit (e.g., 2 becomes 3, and so on), and the last measurement (250) will be dropped from the memory list.

Deleting measurement data

Delete the measurements by following the steps below.

A: Delete one measurement (Example shown below for User 1)

1. Enter the memory mode as described in the previous section "Storing measurement data", find out the measurement you want to delete.
2. Press and hold the Memory button about 3 seconds, the LCD will display "dEL y" and blink.
3. Press the Power/Measurement button  to confirm the deletion, the LCD will display "dEL donE" and then the previous measurement will be displayed.

Note: Before you confirm the deletion, you can press the User or Memory button to switch the LCD display from "dEL y" to "dEL no". Press the Power/Measurement button  to exit delete mode.


B: Delete all measurements (Example shown below for User 1)

1. Enter the memory mode as described in the previous section "Storing measurement data".
2. Press and hold the User and Memory button about 3 seconds, the LCD will display "dEL AL" and blink.
3. Press the Power/Measurement button  to confirm the deletion, the LCD will display "dEL donE".

Note: Before you confirm the deletion, you can press the User and Memory button to switch the LCD display between "dEL AL" and "dEL no". Then press the Power/Measurement button  to exit delete mode.


4. Once deleted, your readings cannot be restored. The LCD will display "----" like the following picture.

Press the Power/Measurement button  to turn off the monitor, otherwise it will power off automatically within 30 seconds.

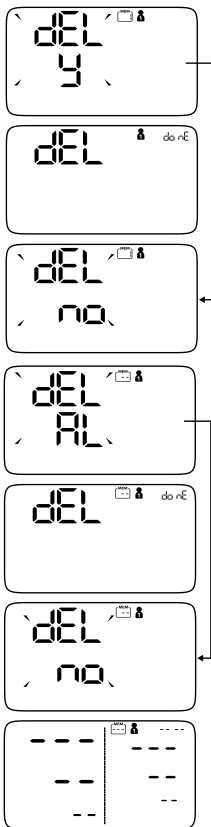Voice Modes and Voice Prompts

The BUA6450 Blood Pressure Monitor includes voice functionality to assist users during use.

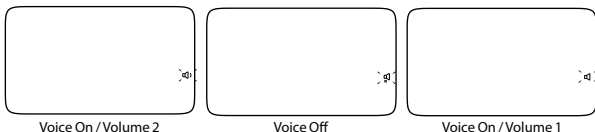Voice On: Enables voice prompts during setup, measurement, and data management.

Voice Off: Disables voice prompts.

Volume Levels: Two volume levels are available when the voice mode is enabled:

Volume 1: Lower volume.

Volume 2: Higher volume.

Voice Prompts During Setup

Setting Date and Time:

When entering setup mode, the monitor will announce: *"Please set up the product before measurement."*

When setting the year: *"Set the year."*

When setting the date format: *"Set the date format. You can set the format as day, month, year or month, day and year."*

When setting the date: *"Set the date."*

When setting the time format: *"Set the time format."*

When setting the time: *"Set the time."*

When setting the voice volume: *"Set the voice volume."*

Voice Settings:

You can enter voice settings by pressing and holding the "USER" button for about 3 seconds.

The monitor will display and announce the current voice mode and volume level.

Voice Prompts During Measurement

Starting Measurement:

If voice mode is enabled, the monitor will announce: *"Starting measurement, please remain still and quiet."*

At beginning of test, voice: *"Select user--user 1/user 2/Guest"*

After measurement, the monitor will broadcast the results, including systolic pressure, diastolic pressure, and pulse rate.

Irregular Pulse or Excessive Body Motion:

If irregular pulse rate or excessive body motion is detected, the monitor will display the respective symbols, and the voice mode will announce the issue.

Voice Prompts During Data Management

Recall Records:

When recalling records, the monitor will announce the average value or the latest record value.

Delete Records:

Delete One Record: The monitor will announce: *"Delete current data."* and *"Successfully deleted."*

Delete All Records: The monitor will announce: *"Delete all data."* and *"Successfully deleted."*

If deletion is canceled, the monitor will announce: *"Cancel."*

Additional Notes

The voice mode can be set during the initial setup or by pressing and holding the "USER" button for about 3 seconds.

The voice prompts are available for User 1 and User 2 but not for the Guest mode.

The voice mode can be turned off if preferred.

This feature is designed to enhance user experience, especially for individuals who may benefit from auditory guidance during operation.

Care and cleaning

In order to get the best performance, please follow the instructions below.

1. Cleaning Process:

Step 1: Make sure to switch off and unplug the device prior to cleaning.

Step 2: Use a soft cloth wetted with soapy water to clean the cuff, and then use a soft cloth wetted with clear water to remove residual soap until there is no visible residual contaminants. Attention shall be paid to avoid liquid invasion into the cuff.

Step 3: Use a dry soft cloth to wipe the cuff, in order to remove residual moisture.

Step 4: Dry the cuff at a well-ventilated place after cleaning.

Step 5: Clean the device with a soft dry cloth.

2. Disinfection Process:

Step 1: Make sure to switch off and unplug the device prior to disinfection.

Step 2: Use a soft cloth wetted with 70% isopropanol to disinfect the cuff for about 3 minutes. Attention shall be paid to avoid liquid invasion into the cuff.

Step 3: Use a clean dry cloth or towel to wipe off the disinfectant until there is no visible residue.

Step 4: Dry the cuff at a well-ventilated place after disinfection.

Suggestion:

Frequency:

For single patient multiple use, it's recommended to clean the device surface once a month or whenever it's necessary.

For multiple patient multiple use, it's recommended to clean the device every time before and after usage. Maintenance procedures shall be taken as per instruction.

The cuff contains a sensitive air-tight inflatable bladder. Handle this carefully and avoid all types of straining through twisting or buckling.

Do not drop the device or treat it roughly in any way. Avoid strong vibrations.

Calibration

This product has been calibrated at the time of manufacture. If used according to the instructions for use, periodic recalibration is not required. If at any time you question the measuring accuracy, please contact our service representative (see last page for contact information).

Errors and troubleshooting

Error	Description	Potential cause and remedy
No power	Display cannot light up	Batteries are depleted. Replace with new batteries. Batteries are inserted incorrectly. Insert the batteries correctly. Adapter is inserted incorrectly. Insert the AC adapter correctly.
High battery	Usb Err shows	USB error. Unplug the adapter to disconnect the device from main power and replace with an authorized adapter.
Low battery	 bAt Lo shows	The battery is too low. Replace with new batteries.
Error message	E 01 shows E 02 or  shows E 03 shows E 04 shows EExx shows  Usb Err shows	The cuff is loose. Re-fit the cuff before measurement. Excessive body motion (such as shaking of the arm with the cuff on or weak Pulse is detected.) Do not move during the measurement and measure again. No pulse is detected. Wear the cuff again before measurement. The measurement failed. Relax for 5 minutes and measure again. Hardware error (XX can be some digital symbol, such as 01, 02, 03, etc.). Turn off monitor and measure again. If EExx still appears on the display, please contact the retailer or our customer service. The voltage of the blood pressure monitor is high. Replace with the authorized adapter.
Warning message	Out shows	Out of measurement range. Relax for a moment and then measure again. If the problem persists, contact your physician.

Warranty

Please read all instructions before attempting to use this device. Please retain the receipt as proof of and date of purchase. The receipt must be presented when making any claim within the relevant warranty period. Any claim under the warranty will not be valid without a proof of purchase.

Your device is warranted for two years (2 years) from date of purchase.

This warranty covers defects in materials or workmanship that occur under normal use; defective devices meeting these criteria will be replaced free of charge.

The warranty DOES NOT cover defects or damage resulting from abuse or failure to follow the user instructions. The warranty becomes void if the device is opened, tampered with, or used with non-Braun branded parts or accessories, or if repairs are undertaken by unauthorized persons.

Accessories and consumables are excluded from any warranty.

For support requests or warranty claims, please contact us by calling 1-800-327-7226 or emailing us at ConsumerRelations@HelenofTroy.com.

Product specifications

Model:	BUA6450
External dimensions:	Approx. 4.26" x 7.14" x 8.12"
Display mode:	Digital LCD V.A. 4.8" x 2.93"
Weight:	Approx. 0.85 lb (Excluding the batteries and cuff)
Measurement mode:	Oscillographic testing mode
Mode of operation:	Continuous operation
Measurement range:	Rated cuff pressure: 0 mmHg~299 mmHg, Measurement pressure: SYS: 60 mmHg ~ 230 mmHg, DIA: 40 mmHg ~ 130 mmHg, Pulse value: (40-199) beat/minute
Accuracy:	Pressure: within ± 3 mmHg, Pulse value: $\pm 5\%$
Operating temperature:	41°F to 104°F
Operating humidity:	15% to 90% relative humidity, non-condensing
Storage/transport temperature:	-04°F to 140°F
Storage/transport humidity:	< 93% relative humidity, non-condensing
Operate atmospheric pressure:	700~1060 hPa
Store atmospheric pressure:	500~1060 hPa
Measurement perimeter of the upper arm:	About 8.6"-17.7", is intended for adult population or those who suffer from diabetes.
Degree of protection:	Type BF applied part
Device Classification:	Battery Powered Mode: Internally Powered ME Equipment, AC Adapter Powered Mode: Class II ME Equipment
Expected Lifetime:	Device: 3 years or 30,000 measurements (may vary based on usage conditions), Cuff: 10000 times, Alkaline battery: About 200-300 times
Types of use/reuse:	Multiple patient multiple use
WARNING:	No modification of this equipment is allowed.
Authorized component:	Please use the authorized adapter (Optional!).
Adapter Type:	Input: 100-240V, 50-60Hz, 0.2A max Output: 1A

If device is not used within specified temperature and humidity ranges the technical accuracy cannot be ensured.

IP21 Protected against solid foreign objects of 12.5 mm diameter and greater.
Protected against vertically falling water drops.

Refer to instruction manual (blue symbol)

Type BF applied part (cuff)

Consult electronic instructions for use

Caution

Operating temperature

Storage temperature

Manufacturer

Date of manufacture

Operating humidity

Storage humidity

Artery mark

Batch code

Serial number

Catalogue number

Medical Device

Importer

Distributor

Unique Device Identifier

Subject to change without notice.

Supply connection: internally powered equipment.

Mode of operation: continuous operation.

The device needs special precautions regarding EMC. Portable and mobile RF communications equipment can affect medical electrical equipment. **WARNING:** Use of this equipment adjacent to other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally. **WARNING:** Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30cm to any part of the device. Otherwise, degradation of the performance of this equipment could result. **WARNING:** use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation. For more information, see consumer service numbers or website.

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